

Angio cpet - with electrical adjustable wall fixation

Modern ergometer with multifunctional applications



Highlights

High standards

Lode is a socially and environmentally responsible company. All Lode products are RoHS/WEE compliant and Lode is ISO 9001:2003, ISO and 13485:2008 certified. All medical products comply to MDD 93/42/EEC, incl. IEC 60601-1.

Easy to operate

For Lode products this means:

- easy to connect
- easy to move around
- easy user interface

Reliable and reproducible stress tests

The experience of professionals who calibrate many ergometers shows that the Lode ergometers are the most reliable across the complete workload and rpm range and still within specifications even after many years of intensive use.

Additional features with PCU

Besides the possibility to program 24 protocols easily, this control unit offers the following features:

- better monitoring because of the additional and larger display
- a perfect combination with BPM
- possibility to measure SpO2



Angio cpet - with electrical adjustable wall fixation



Modern ergometer with multifunctional applications

The wall fastening consists of a vertical rail with hang-up mechanism which makes it possible to move the Angio up and down electrically over a range of 70 cm. Handgrips are standard included in this combination. The Angio is an ergometric unit that can be used for both arm and supine ergometry. Its compact design makes it universally applicable for ergometric studies in those sectors in which standard ergometry cannot be used. The Angio operates independent of pedaling speed in the range of 7 - 1000 watt. The Angio cpet is standard supplied with a communication module and can therefore be easily controlled by all known stress ECG and pulmonary devices in the world. The workload, rpm and time can be readout from the 3,5" colour display.

Features



Accurate over a long period of time

The Lode ergometers are supplied with an electro-magnetic braking mechanism of Lanooy (eddy current). The biggest advantage of this braking system compared to a friction braking system is the absolute accuracy and the accuracy over time. Moreover, friction braking systems have more wearing parts.

**7
watt**

Extreme low start-up load

The extreme low start-up load of 7 watts and the adjustability in small steps of 1 watt make this ergometer perfectly suitable for many different applications. The standard control unit shows multiple ergometry parameters and you can determine your specific default setting and start-up menu.

**1
watt**


Small adjustment steps

The workload of the Lode ergometers is adjustable in steps of only 1 watt. Depending on your wishes, the test operator or the test subject can adjust the workload. The steps of 1 watt are possible in the manual mode as well as within protocols.



RS232 connectivity

RS232 ports enable connectivity to most ECG and ergospirometry devices as well as PC's.

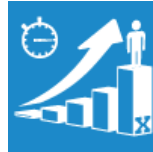

LEM compatible

This product can be used with Lode Ergometry Manager (LEM) software to manage data and to apply specific protocols when a Communication card or Network card is present



LCRM compatible

This product can be used with Lode Cardiac Rehabilitation Manager software (LCRM)

Lode













FOR LIFE | UNDERSTANDING
MOVEMENT & PERFORMANCE

Angio cpet - with electrical adjustable wall fixation



Modern ergometer with multifunctional applications

Angio cpet - with electrical adjustable wall fixation can a.o be extended with the following options:

<p>Pedal shoes (pair)</p> <p>Extra stability during cycling</p>  <p>Partnumber: 917803</p>	<p>Pedal shoes pediatric (pair)</p> <p>Pedal shoes for children</p>  <p>Partnumber: 917833</p>	<p>Pedal shoes extra large (pair)</p> <p>For large feet sizes</p>  <p>Partnumber: 917834</p>	<p>Adjustable cranks</p> <p>Optimal force application</p>  <p>Partnumber: 928804</p>	<p>USB to Serial converter</p> <p>Easy connection</p>  <p>Partnumber: 226012</p>
<p>Handgrip fixation set</p> <p>Handgrip fixation set</p>  <p>Partnumber: 917822</p>	<p>Handgrip fixation set pediatric</p> <p>Versatile ergometry</p>  <p>Partnumber: 917823</p>	<p>RS232 cable</p> <p>Easy connection</p>  <p>Partnumber: 930911</p>	<p>Control Unit with 7" touch screen for ergometer</p> <p>Multifunctionality</p>  <p>Partnumber: 945834</p>	<p>Programmable Control Unit with 7" Touchscreen for ergometer</p> <p>Programmable</p>  <p>Partnumber: 945835</p>
<p>SpO2 for control unit with touch panel (bicycle)</p> <p>Saturation and heart rate</p>  <p>Partnumber: 945823</p>	<p>Heart rate large distance</p> <p>Heart rate in beats per minute</p>  <p>Partnumber: 945830</p>			

Angio cpet - with electrical adjustable wall fixation



Modern ergometer with multifunctional applications

Specifications

Workload

Minimum load	7 W
Maximum peak load	1000 W
Minimum load increments	1 W
Maximum continuous load	750 W
Hyperbolic workload control	✓
Linear workload control	✓
Fixed torque workload control	✓
Maximum rpm independent constant load	150 rpm
Minimum rpm independent constant load	30 rpm
Optional heart rate controlled workload	✓
Electromagnetic "eddy current" braking system	✓
Dynamic calibration	✓

Accuracy

Workload accuracy from 7 to 100 W	3 W
Workload accuracy from 100 to 500 W	3 %
Workload accuracy from 500 to 1000 W	5 %

Comfort

Adjustability range height	700 mm
----------------------------	--------

User Interface

Readout Distance	✓
Readout RPM	✓
Readout Heartrate	✓
Readout target HR	✓
Readout Energy	✓
Readout Torque	✓
Readout Time	✓
Readout Power	✓
Set Display	✓
Set Resistance	✓
Set P-Slope	✓
Set Mode	✓
Manual operation mode	✓
Preset protocol operation mode	✓
Terminal operation mode	✓
External control unit	✓
Selfdesigned protocol operation mode	✓

Dimensions

Product length (cm)	68 cm	26.8 inch
Product width (cm)	42 cm	16.5 inch
Product height	104 cm	40.9 inch
Product weight	55 kg	121.3 lbs

Power requirements

V AC	100 - 240 V
Phases	1
Frequency	50/60 Hz
Power consumption	250 W
Power cord IEC 60320 C13 with CEE 7/7 plug	✓
Power cord NEMA	✗

Standards & Safety

IEC 60601-1:2012	✓
ISO 13485:2016 compliant	✓
ISO 9001:2015 compliant	✓

Certification

CE class Im according to MDD93/42/EEC	✓
CE class of product with optional SpO2	Ila
CE class of product with optional BPM	Ila
CB according to IECCE CB - pending	✓

Order info

Partnumber: 967903

*Specifications are subject to change without notice.