

# Angio rehab - with Fixation Set for Instrument Rail

Modern ergometer with multifunctional applications



## Highlights

### High standards

Lode is a socially and environmentally responsible company. All Lode products are RoHS/WEE compliant and Lode is ISO 9001:2015, and ISO 13485:2016 certified. All medical products comply to MDD 93/42/EEC, incl. IEC 60601-1.

### Reliable and reproducible stress tests

The experience of professionals who calibrate many ergometers showed that the Lode ergometers are the most reliable across the complete workload and rpm range and still within specifications even after many years of intensive use.

### Various test modes

Besides the hyperbolic (rpm-independent) mode that is used most of the time, the standard control unit offers several other test modes, like the fixed torque mode and the linear mode. These modes can be used in both manual and terminal mode.

### Multifunctional

The ergometer can be used in various ergometry settings, enabling a multifunctional deployment.

### Easy to clean

The design of the ergometer and the housing material make it very easy to clean.



# Angio rehab - with Fixation Set for Instrument Rail



## Modern ergometer with multifunctional applications

In this system, the Angio is supplied with a multifunctional fixation set. With this set, the Angio ergometer is equipped with pedal shoes and can be mounted easily on all Imaging / Physiotherapy / Examination tables that have a standard instrument side rail of 25x10 mm. The adjustability range is 625 - 903 mm (width) and 20-120 mm (height). \*

The Angio rehab is standard supplied with a 7" control unit with touch screen. Thanks to the built-in network module, the ergometer can be connected to the Lode Cardiac Rehab Manager Software, Lode Rehab Manager or the Lode Ergometry Manager.

A USB A-B cable for service purposes only will be standard delivered with the product. To connect LEM or LCRM you need a special interface cable that can be ordered under part number 930930.

Note: always check the width of the examination table on which the ergometer is going to be installed. There is a fixation set available for tables that are smaller (#967820)

# Angio rehab - with Fixation Set for Instrument Rail



Modern ergometer with multifunctional applications

## Features

**7  
watt**

### Extreme low start up load

The extreme low start-up load of 7 watts and the adjustability in small steps of 1 watt make this ergometer perfectly suitable for many different applications. The standard control unit shows multiple ergometry parameters and you can determine your specific default setting and start-up menu.



### Accurate over a long period of time

The Lode ergometers are supplied with an electro-magnetic braking mechanism of Lanooy (eddy current). The biggest advantage of this braking system compared to a friction braking system is the absolute accuracy and the accuracy over time. Moreover, friction braking systems have more wearing parts.

**1  
watt**


### Small adjustment steps

The workload of the Lode ergometers is adjustable in steps of only 1 watt. Depending on your wishes, the test operator or the test subject can adjust the workload. The steps of 1 watt are possible in the manual mode as well as within protocols.



### RS232 connectivity

RS232 ports enable connectivity to most ECG and ergospirometry devices as well as PC's.

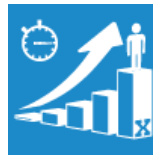

### LEM compatible

This product can be used with Lode Ergometry Manager (LEM) software to manage data and to apply specific protocols when a Communication card is present



### LCRM compatible

This product can be used with Lode Cardiac Rehabilitation Manager software (LCRM)



### Versatile controls

### Additional features with PCU

Besides the possibility to program 24 protocols easily, this control unit offers the following features:

- better monitoring because of the additional and larger display
- a perfect combination with BPM
- possibility to measure SpO2



### Versatile Interfacing

Various interface protocols guarantee perfect communication with all commonly known stress ECG and spirometry equipment

# Angio rehab - with Fixation Set for Instrument Rail



Modern ergometer with multifunctional applications



## Interconnectivity between Lode products

Connecting Lode products has never been easier! Lode rehab and sports products have a standard Network card:

- To be able to connect the first product to the PC with L(C)RM a Lode proprietary network to PC cable is needed (#930930). This cable is standard included with Lode Rehab Software.
- From the second product onward products can be connected to the previous one, creating a bus network configuration;
- The last product always needs a termination plug to avoid interference and loss of data. Therefore all products with such a network card come with a termination plug.

### Benefits




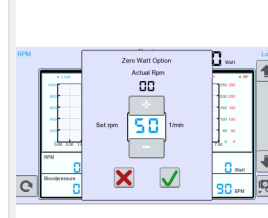
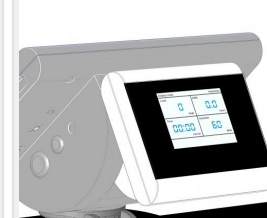









- Lossless data connection
- High bandwidth
- No interference of COM ports
- Daisy chain connection
- Full access of all data in the product to LCRM

# Angio rehab - with Fixation Set for Instrument Rail



Modern ergometer with multifunctional applications

Angio rehab - with Fixation Set for Instrument Rail can a.o be extended with the following options:
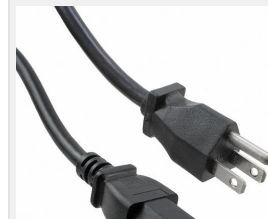
<p>Pedal shoes (pair)</p> <p>Extra stability during cycling</p>  <p>Partnumber: 917803</p>	<p>Pedal shoes pediatric (pair)</p> <p>Pedal shoes for children</p>  <p>Partnumber: 917833</p>	<p>SpO2 for control unit with touch panel (extra long cable)</p> <p>Oxygen saturation</p>  <p>Partnumber: 945822</p>	<p>0-Watt start-up system</p> <p>Lowest possible startup power</p>  <p>Partnumber: 960805</p>	<p>Colour Display 3.5" - 2nd screen</p> <p>Multifunctionality</p>  <p>Partnumber: 945819</p>
<p>Network module connection cable</p> <p>Easy connection</p>  <p>Partnumber: 930930</p>	<p>SpO2 for control unit with touch panel - ordered afterwards</p> <p>Ordered afterwards</p>  <p>Partnumber: P945822</p>	<p>Add program function to 7" touch screen for ergometer</p> <p>Easily programmable</p>  <p>Partnumber: U945835</p>	<p>Blood Pressure with ECG trigger for bicycle ergometer with ECG trigger</p>  <p>Partnumber: 945828</p>	<p>Wheeled frame for Angio Ergometer w. fixation set</p> <p>Easy installation of the Angio ergometer</p>  <p>Partnumber: 967800</p>
<p>Small Fixation Set for Instrument Rail</p> <p>For smaller examination tables of</p>  <p>Partnumber: 967820</p>	<p>Bluetooth Smart heart rate</p> <p>Heartrate available within an extreme wide</p>  <p>Partnumber: 945833</p>	<p>Adjustable cranks</p> <p>Optimal force application</p>  <p>Partnumber: 960804</p>	<p>Custom pack Powercable 100 - 127 V for ergometers</p> <p>Easy connection</p>  <p>Partnumber: 945892</p>	

# Angio rehab - with Fixation Set for Instrument Rail



Modern ergometer with multifunctional applications

## Specifications

### Workload

Workload range fixed torque	0,1 - 70 Nm
Minimum load	7 W
Maximum peak load	1000 W
Minimum load increments	1 W
Maximum continuous load	750 W
Hyperbolic workload control	✓
Linear workload control	✓
Fixed torque workload control	✓
Maximum rpm independent constant load	150 rpm
Minimum rpm independent constant load	30 rpm
Optional heart rate controlled workload	✓
Electromagnetic "eddy current" braking system	✓
Dynamic calibration	✓

### Accuracy

Workload accuracy from 7 to 100 W	3 W
Workload accuracy from 100 to 500 W	3 %
Workload accuracy from 500 to 1000 W	5 %

### Comfort

Adjustability range height	20 - 120 mm	
Adjustability range width	625 - 903 mm	24.6 kJ
Pedal shoes	✓	

### User Interface

English user interface	✓
Chinese user interface	✓
Croatian user interface	✓
Czech user interface	✓
Danish user interface	✓
Dutch user interface	✓
Finnish user interface	✓
French user interface	✓
German user interface	✓
Greek user interface	✓
Hungarian user interface	✓
Italian user interface	✓
Japanese user interface	✓
Korean user interface	✓
Latvian user interface	✓
Lithuanian user interface	✓
Norwegian user interface	✓
Polish user interface	✓
Portugese user interface	✓
Romanian user interface	✓
Russian user interface	✓
Spanish user interface	✓
Swedish user interface	✓
Turkish user interface	✓
Ukrainian user interface	✓
Readout Distance	✓
Readout RPM	✓
Readout Heartrate	✓
Readout target HR	✓
Readout Energy	✓
Readout Torque	✓
Readout Time	✓
Readout Power	✓
Set Display	✓
Set Resistance	✓
Set P-Slope	✓
Set Mode	✓
Manual operation mode	✓
Preset protocol operation mode	✓
Analog operation mode	✓
Terminal operation mode	✓
External control unit	✓
Selfdesigned protocol operation mode	✓

# Angio rehab - with Fixation Set for Instrument Rail



Modern ergometer with multifunctional applications

## Dimensions

Maximum table width	920 mm	36.2 inch
Minimum table width	625 mm	24.6 inch
Product length (cm)	54 cm	21.3 inch
Product width (cm)	68 cm	26.8 inch
Product height	73 cm	28.7 inch
Product weight	40 kg	88.2 lbs

## Power requirements

V AC	100 - 240 V
Phases	1
Frequency	50/60 Hz
Power consumption	160 W
Power cord IEC 60320 C13 with CEE 7/7 plug	✓
Power cord NEMA	✗

## Standards & Safety

IEC 60601-1:2012	✓
ISO 13485:2016 compliant	✓
ISO 9001:2015 compliant	✓

## Certification

CE class Im according to MDD93/42/EEC	✓
CE class of product with optional SpO2	Ila
CE class of product with optional BPM	Ila
CB according to IEC60957 CB	✓

## Order info

Partnumber: 967925

*\*Specifications are subject to change without notice.*